

CAGENIX PRODUCT GUIDE



C A G E N I X™

Beyond the Bar: 14,000 arches and still growing strong.

Cagenix set out to manufacture the finest frameworks in the world. Now, we're innovating across the full-arch spectrum.

With AccuFrame, Cagenix perfected the fixed-hybrid using sophisticated CAD/CAM technology and a patented proprietary imaging process.

Then, innovations like AccuFrame IC and Cagenix 360 advanced what's possible in restorative technology, establishing new benchmarks for quality, versatility, and durability in completed restorations.

Now, with patents pending for our new BarZero conversion and provisional systems, Cagenix is poised to unveil yet another in a history of ongoing advances from the first name in frames.



Introducing BarZero

The NEW frameless conversion and provisional system

What happens when the company that revolutionized the framework sets their sights on improving the world of conversions and provisionals? Conversions get faster and more reliable. Your provisional-to-final process is streamlined. Patient satisfaction improves along with your bottom line.

With two options available, BarZero Express provides a bulletproof day-of-surgery provisional that virtually eliminates repair visits, and BarZero LTP satisfies your need for high-quality long-term provisional.

AccuFrame: The original and still the best

We created the first colorized titanium fixed-hybrid framework, and then we perfected it.

It's been years since Cagenix took the fixed-hybrid world by storm with AccuFrame & AccuFrame Plus setting the industry standard for quality and accuracy. Today, you can still experience the AccuFrame difference with our custom-milled titanium frameworks.



Our premium restorative systems: a cut above the rest

After perfecting the fixed-hybrid, we set our sights on advancing restorative science through innovation.



AccuFrame IC

Ceramic individual crowns on our titanium frame
AccuFrame IC restorations utilize individual crowns in a choice of ceramic materials, including high-strength zirconia, e.Max® or even layered porcelain, on our award-winning frameworks.



Cagenix 360


Our unique overlay-on-universal-frame system
An overlay-on-frame system built around a universal Titanium bar designed to support a range of interchangeable restorations for the ultimate in customization and reparability.




One frame. Endless possibilities.

The game-changing overlay-on-frame restorative system.

Cagenix 360 features a custom milled overlay structure that mates to a milled Titanium framework. The unique framework is modular in design, allowing easy replacement and interchangeability within the current and future overlay options in the system. The overlays are milled from a range of materials including PMMA and Zirconia, offering flexibility both upon delivery and for future upgrade considerations.

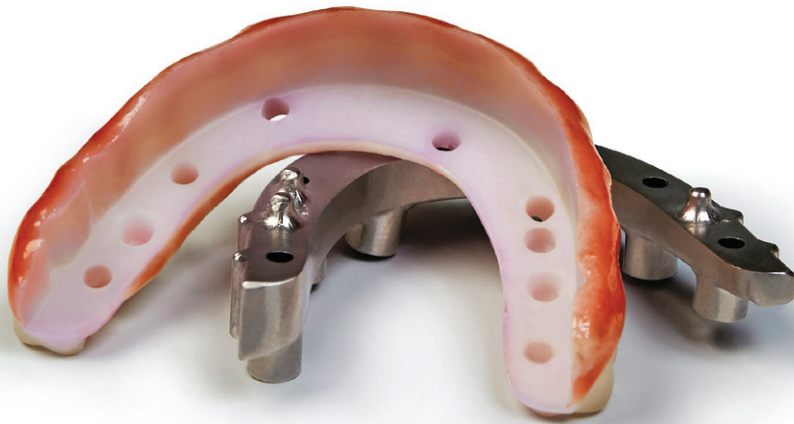

- As little as 8mm restorative space required
- Lifetime framework warranty
- Easy repair/replacement in as little as 1 visit
- Choose from full wrap or metal-base options
- Supports any current or future millable material
- Supports industry's widest range of implant platforms
- Receive restoration in as little as 14 business days

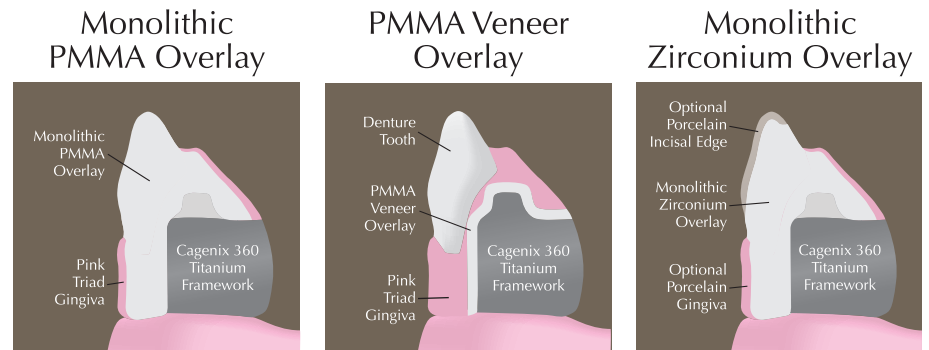

Select from 3 overlay options to customize the restoration to your needs

ACCUFRAME IC™

INDIVIDUAL CROWN SYSTEM



Ceramic individual crowns on our titanium framework

AccuFrame IC restorations utilize individual crowns in a choice of ceramic materials, including high-strength zirconia, e.max® or even layered porcelain, on our award-winning AccuFrame or AccuFrame Plus titanium frameworks. The result? Guaranteed fit, superior appearance and long-lasting performance.

- Stronger than traditional fixed-hybrid restorations: Virtually eliminates tooth wear and breakage.
- Easy to get started: No investment in lab equipment
- Simple processing: Requires only your order and set-up
- Supports any current or future millable material
- Receive restoration in as little as 14 business days

ACCUFRAME®



The gold-standard for fixed-hybrid milled frameworks.

With the AccuFrame and AccuFrame Plus frameworks, Cagenix has set the standard for quality and service in the world of fixed-hybrid restorations.

- Lifetime warranty for all titanium restorations
- Guaranteed fit, first-time, every-time
- Aesthetic “Plus” pink-frame treatment eliminates need for opaquing
- Support for nearly all implant connections
- Receive restoration in as little as 10 business days

BARZERO

Cagenix: Going beyond the bar

Introducing the barless conversion and provisional system
Whether it's an instant conversion or long-term provisionalization, patent-pending BarZero restorations offer improved durability, aesthetics, and convenience over competing options.



BARZERO EXPRESS

Milled PMMA Conversion Prosthesis
The faster, more durable immediate conversion prosthesis for your patients

- Install on day-of-surgery in about an hour
- “Trough” design picks up copings intra-orally then secures them with acrylic
- Improved strength over traditional options
- Virtually eliminate reduce repair incidents
- Available with a palate and flanges if desired



BARZERO LTP

Milled PMMA Long-Term Provisional
Durable monolithic PMMA barless restoration provides excellent value

- Ideal interim prosthesis to temporize the patient during osseointegration
- Arrives complete and ready to install
- Improved strength over traditional options
- Monolithic design reduces repair incidents
- Improved aesthetics & patient satisfaction

Supported Connections

We know you have enough to worry about without adding the concern of platform support to the mix. That is why we strive to support one of the widest ranges of implant and abutment systems in the industry.

If you don't see a connection listed here, give us a call. We're adding new connections all the time, and support many legacy systems.

(866) 964-5736

AB Dental	Angular adapter	MegaGen	Anyridge Fixture (3.5/4.0/6.0) ExFeel External Hex Fixture (5.0) EZ Plus Analog (3.3/4.0) Anyone Fixture (Small/Regular) Multi-Unit Abutment
BioHorizons	Abutment for Screw (4.0/5.0) External Hex (3.5/4.0/5.0) Tapered Internal Implant (3.0/3.5/4.5/5.0) Multi-Unit Abutment	MIS	TransOcta Abutment Internal Hex (NP/RP/WP) Old Multi-Unit Abutment New Multi-Unit Abutment
Biomet 3i	Conical Abutment (3.4/4.1/5.0/6.0) Certain Fixture (3.4/4.1/5.0) External Hex (3.4/4.1/5.0/6.0) Standard Abutment Low Profile Abutment IOL Abutment	Neoss	Fixtures Access Abutments
Camlog	Internal Fixture (3.3/3.8/4.3/5.0) Bar Abutment (4.3/6.0) Conelog Fix (3.3/3.8/4.3/5.0) Vario SR Abutment (4.3 only)	Nobel Biocare	Multi-Unit Abutment Branemark External Hex Fixture (NP/RP/WP) Replace Internal Fixture (NP/RP/WP/6.0) Active Fixture (3.0/NP/RP) Replace Conical Connection (NP/RP) Branemark WP Multi-Unit Abutments Sterioss PME Abutment Novum External Hex Standard Abutment
Dentis	S-Clean Fixture Level Multi-Unit Abutment	Southern	External Hex Fixture (3.25/3.75/4.0/5.0/) Internal Tri-Nex (3.5/4.3/5.0/6.0) Compact Conical Abutment (Narrow and Wide) Octabutment for 4.8 Internal Octagon Fixture
Dentsply / Astra	UniAbutment (20°/45°) OsseoSpeed Fixtures (3.0, 3.5/4.0, 4.5/5.0) EV Fixture (3.0/3.6/4.2/4.8/5.4) EV Unit Abutment (Straight Only) Ankylos Implant Fixture Balance Base Abutment (Narrow and Standard) Frialit Fixture (3.3/3.8/4.5/5.5)	Straumann	Synocta 1.5mm Abutment (RN/WN) Soft Tissue (NN/RN/WN) Bone Level Fixture (NC/RC/NNC 3.5) NC/RC Screw Retained Abutment (SRA) Multi-Base Abutment (NC/RC/BC)
Glidewell	Inclusive Fixture (3.5/4.5)	Sybron	Octa Abutment System Pitt-Easy Fixture Level (3.75)
Hiossen	HS Octa Abutment Multi-Angled Esthetic-low Abutment Covertable Abutment (4.0/4.8/6.0) Fixture (Mini/Regular)	Thommen	SPI VARIO Multi-Abutment Fixture (3.5/4.0/4.5/5.0/6.0)
Hi-Tec	LGI Plus Fixture Logic 3.5/4.3/5.0 Fixture	Zimmer	Eztetic (3.1) Tapered Screw Vent (3.5/4.5/5.7) Tapered Abutment Advent Fixture Spectra Cone Abutment Calcitek Integral Omilok Paragon Fixture (3.5) Shouldered Abutment for Omnilok TSI Abutment (TSI/TSI1) Spline Fixture (4.0/5.0)
Implant Direct / Attachments Intl.	Overdenture Abutment Legacy Fixtures Tri-Lobe Fixtures Screwplant Direct (3.7/4.7/5.7) Screw Indirect Implant UMA Abutment (Attachments Int) Multi-Use Abutment (attachments Int)		
Keystone	Fixed Detachable Abutment Indexed Abutment (RDS/WDS) Prima Implant (SD 3.5/RD 4.1/ WDS 5.0) Multi-Unit Abutment (RD/WD) Renova (SDI) Restore (RD/SD/WD) OCT Fixture (4.8/6.5)		

Submitting your case

What to send, where to send it, and what to expect.

Below is a good outline of what we require when submitting a case, as well as the typical turn-around times for each of our products. We will generally send a 3D design approval file via e-mail within 3 days of our receiving the case. If you would like a call from one of our engineers prior to designing your case, just let us know. All necessary documents are available for download from <http://cagenix.com>.

Cagenix Products	Completed Order Form	Clean Verified Master Cast	Diagnostic Wax-Up	Screw-Retained Diagnostic Wax Up	Articulator w/ Opposing Arch	Verification Jig	Screws	Turnaround*
AccuFrame	Y	Y	Y	N	N	O	N	10 Business Days
AccuFrame IC	Y	Y	N	Y	Y	O	N	14 Business Days
Cagenix 360	Y	Y	N	Y	Y	O	N	14 Business Days
BarZero	Y	Y	Y	N	N	O	N	10 Business Days

*Delays in design approval, or requests for lab work like gingival or crown processing may require additional time
Cases submitted without required components may be subject to delays or additional labwork charges.

Y = Required **O = Optional** **N = Do Not Include**

Submit cases to: Cagenix Inc. • 1610 Century Center Pkwy • Bldg. F, Suite 102 • Memphis, TN 38134
Product questions or other inquiries: Call us at (866) 964-5736 or visit <http://cagenix.com>

